

MEDICAL

Medical Vacuum

Models | A3052MV to A3303MV

Designed for use in critical medical applications to remove solid, liquid and bacterial contamination in vacuum systems, Walker Filtration's Medical Vacuum Filters guarantee a safe and reliable product that is trusted by hospitals worldwide.

Situated on the suction side of a vacuum pump, the New Alpha Series of Medical Vacuum Filters are essential to avoiding damage to vacuum pumps and preventing potentially hazardous biological contagions from being exhausted into the surrounding environment.



NEW Filtration Technology Filtration efficiency in excess of 99.9999%



Product Safety in Mind
Easily removable sterilizable drain
flask and differential pressure
monitor supplied as standard



Push Fit Element Design
Quick and easy maintenance
with unique push fit
element design

- High Efficiency Medical Grade Filter Element Custom engineered filter media and deep pleat element technology provides minimal pressure loss and filtration efficiencies in excess of international medical gas standards
- International Medical Gas Standards Fully compliant with global medical gas pipeline systems standards ISO 7396-1, NFPA 99, HTM 02-01 and AS 2896
- Quick and Easy Maintenance Unique push fit element design allows for easy maintenance, significantly reducing service time and contact time for service engineers with contaminated filter elements
- Advanced Filtration Technology Low pressure loss borosilicate glass microfiber-media and open cell
 reticulated foam pre-filtration layer captures particulates, bacteria and liquid aerosols reducing energy
 consumption and overall system costs for low total cost of ownership
- Corrosion Protection Internal and external electrophoretic painting followed by a tough exterior polyester powder coating
- Product Safety in Mind Quality design and build. Guaranteed safe housing closure with rotational safety stop







For further information please visit www.walkerfiltration.com



Third party tested and validated in accordance with NFPA 99, HTM 02-01, ISO 7396-1 and AS 2896









Technical Specification

Filter model	Pipe size inches	Maximum Rated Flow at an Operating Vacuum (Suction) Pressure of 19"Hg (475 mmHg) [(285 torr)]						Dimensions inches (mm)				Weight	Weight	Element
		Free Air Aspirated (FAA)			Rarified Air Flow							lbs	Kg	model
		SCFM	NL/min	Nm³/hr	CFM	L/min	m³/hr	Α	В	С	D			
A3052MV	1/2	6.7	190	11.4	17.7	502	30.1	5.00 (127)	1.26 (32)	11.22 (285)	3.15 (80)	3.7	1.7	E30612MV
A3071MV	3/4	8.3	235	14.1	21.9	621	37.3	5.00 (127)	1.26 (32)	11.22 (285)	3.15 (80)	3.7	1.7	E30612MV
A3101MV	1	9.7	275	16.5	25.7	727	43.6	5.00 (127)	1.26 (32)	11.22 (285)	3.15 (80)	3.7	1.7	E30612MV
A3102MV	1	17.1	485	29.1	45.3	1282	76.9	5.00 (127)	1.26 (32)	14.61 (371)	3.15 (80)	4.4	2.0	E30621MV
A3122MV	11/4	33.5	950	57.0	88.7	2512	150.7	6.69 (170)	2.09 (53)	20.00 (508)	3.94 (100)	10.8	4.9	E30831MV
A3151MV	11/2	37.1	1050	63.0	98.0	2776	166.6	6.69 (170)	2.09 (53)	20.00 (508)	3.94 (100)	10.8	4.9	E30831MV
A3201MV	2	38.8	1100	66.0	102.7	2908	174.5	6.69 (170)	2.09 (53)	20.00 (508)	3.94 (100)	10.8	4.9	E30831MV
A3202MV	2	54.7	1550	93.0	144.7	4098	245.9	6.69 (170)	2.09 (53)	27.87 (708)	3.94 (100)	12.1	5.5	E30850MV
A3251MV	21/2	92.9	2630	157.8	245.6	6953	417.2	8.66 (220)	2.76 (70)	28.98 (736)	3.94 (100)	23.1	10.5	E31140MV
A3301MV	3	105.9	3000	180.0	280.1	7932	475.9	8.66 (220)	2.76 (70)	28.98 (736)	3.94 (100)	23.1	10.5	E31140MV
A3302MV	3	118.3	3350	201.0	312.8	8857	531.4	8.66 (220)	2.76 (70)	33.74 (857)	3.94 (100)	25.4	11.5	E31160MV
A3303MV	3	125.4	3550	213.0	331.4	9386	563.1	8.66 (220)	2.76 (70)	39.57 (1005)	3.94 (100)	27.6	12.5	E31175MV

NOTE: Rated flows are stated at NFPA 99 conditions. For flows at other operating vacuum pressures, please consult Walker Filtration Inc. Operating vacuum (suction) pressure is stated at the filter outlet (i.e. vacuum pump or suction side) Standard (reference) atmosphere: 101.325 kPa (1013.25 mbar(a)), 68°F (20°C)




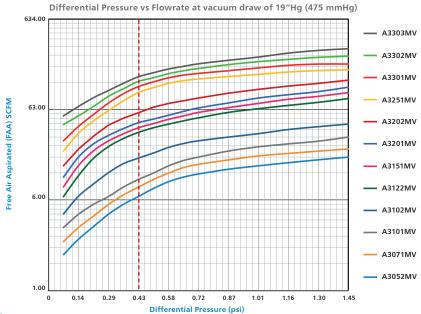

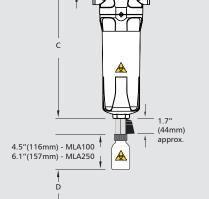

(71mm)

A3052MV - A3303MV

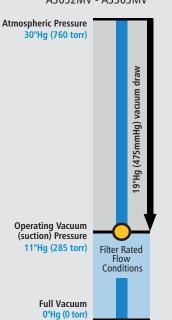

Technical notes

- Direction of flow is outside to inside through the element.
- Filter elements should be replaced at least every 6 months.
- Differential pressure gauges (65DPG250B) are fitted to all models as standard. Volt free contact options are available see price guide.
- Manual drain valves (MDVE25) are fitted to all models. Sterilizable glass drain flasks are supplied as standard, 100ml for models A3052MV to A3101MV and 250ml for models A3102MV to A3303MV.
- Threaded connections are NPT to ANSI/ASME B1.20.1.









